



 **i-gel**<sup>®</sup>

The supraglottic airway with a non-inflatable cuff



**Airway Management** ▪ Airway Devices



## Airway management has evolved

Introducing the i-gel<sup>®</sup>: a revolutionary single use supraglottic airway from Intersurgical.



### i-gel<sup>®</sup> and natural airway management

The i-gel<sup>®</sup> is a truly unique single use, latex and PVC free airway device, representing the culmination of years of extensive research and development. Everything about the i-gel<sup>®</sup> has been designed to work in perfect unison with the anatomy; the i-gel<sup>®</sup> design was inspired by the physiology of the perilaryngeal framework itself – airway management as nature might have intended.

### i-gel<sup>®</sup> mirrors the anatomy

The shape, softness and contours accurately mirror the perilaryngeal anatomy to create the perfect fit. This innovative concept means that no cuff inflation is required. The i-gel<sup>®</sup> works in harmony with the patient's anatomy so that compression and displacement trauma are significantly reduced or eliminated.



### The non-inflatable cuff

i-gel<sup>®</sup> gets its name from the soft gel-like material from which it is made. It is the innovative application of this material that has enabled the development of a unique non-inflatable cuff. This key feature means insertion of i-gel<sup>®</sup> is easy, rapid and consistently reliable.

### The simple, safe and rapid solution

i-gel<sup>®</sup> is incredibly easy to use. A proficient user can achieve insertion of the i-gel<sup>®</sup> in less than 5 seconds. With no inflatable cuff, i-gel<sup>®</sup> provides a safe and rapid airway management solution.




### Adults

Adult i-gel<sup>®</sup> is indicated for securing and maintaining a patient airway in routine and emergency anaesthetics of fasted patients, during spontaneous or intermittent positive pressure ventilation (IPPV), during resuscitation of the unconscious patient, and as a conduit for intubation under fiberoptic guidance in a known or unexpectedly difficult intubation, by personnel who are suitably trained and experienced in the use of airway management techniques and devices.

### Accurate and natural positioning

The i-gel<sup>®</sup> accurately and naturally positions itself over the laryngeal framework, providing a reliable perilaryngeal seal without the need for an inflatable cuff.

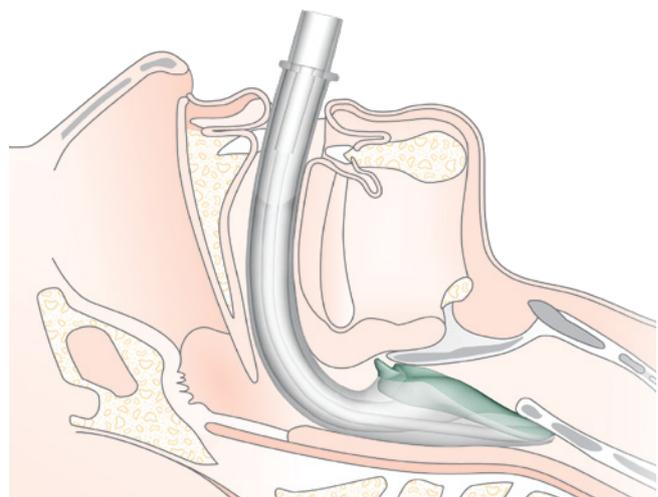

### Paediatrics

i-gel<sup>®</sup> is now available in four paediatric as well as three adult sizes, making it applicable for use with patients between 2–90+kg. Paediatric i-gel<sup>®</sup> is indicated for securing and maintaining a patient airway in routine and emergency anaesthetics for operations of fasted patients during spontaneous or intermittent positive pressure ventilation (IPPV).

### Additional information available

An i-gel<sup>®</sup> user guide, clinical study material and other support documentation is available to download from the i-gel<sup>®</sup> website at [www.i-gel.com](http://www.i-gel.com).

 Video available at [www.intersurgical.com](http://www.intersurgical.com)



## Features and benefits

The i-gel<sup>®</sup> has a host of features that provide significant benefits to the patient and the clinician.

### 15mm connector

Reliable connection to any standard catheter mount or connection

### Proximal end of gastric channel

### Clearly displayed product information

For quick easy reference. Includes confirmation of size and weight guidance

### Position guide (adult sizes only)

Easy confirmation of optimum insertion depth



### Gastric channel

The i-gel<sup>®</sup> incorporates a gastric channel (except size 1). It provides an early warning of regurgitation, allows for the passing of a nasogastric tube to empty the stomach contents and facilitates venting



### Integral bite block

Reduces the possibility of airway channel occlusion

### Buccal cavity stabiliser

Aids insertion and eliminates the potential for rotation


### Epiglottic rest

Reduces the possibility of epiglottic 'down folding' and airway obstruction

### The non-inflatable cuff

Made from a unique soft gel-like material allowing ease of insertion and reduced trauma

### Distal end of gastric channel

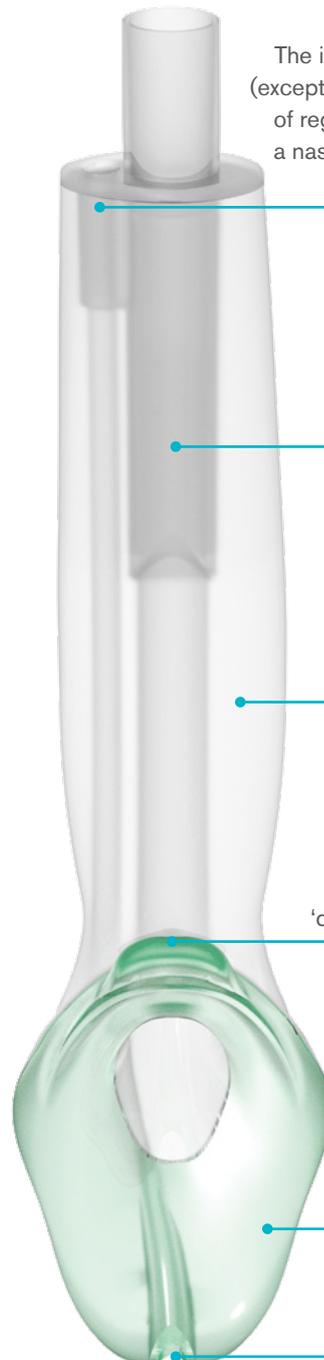

The adult sizes of i-gel<sup>®</sup> can be used as a conduit for intubation under fiberoptic guidance in a known or unexpectedly difficult intubation

### Innovative packaging

The i-gel® supraglottic airway is supplied in a fully recyclable protective cradle or cage pack. This unique packaging protects the i-gel® in transit and ensures that it maintains its anatomical shape. i-gel® is available in seven sizes.



Code	Description	Size	Weight	Box Qty.
8205000	i-gel®, supraglottic airway	5 Large adult	90+kg	25
8204000	i-gel®, supraglottic airway	4 Medium adult	50–90kg	25
8203000	i-gel®, supraglottic airway	3 Small adult	30–60kg	25
8225000	i-gel®, supraglottic airway	2.5 Large paediatric	25–35kg	10
8202000	i-gel®, supraglottic airway	2 Small paediatric	10–25kg	10
8215000	i-gel®, supraglottic airway	1.5 Infant	5–12kg	10
8201000	i-gel®, supraglottic airway	1 Neonate	2–5kg	10

Make an enquiry

Visit the i-gel® website [www.i-gel.com](http://www.i-gel.com)

Sterile

#### References

1. CD Deakin, JP Nolan, J Soar, K Sunde, RW Koster, GB Smith, GD Perkins. European Resuscitation Council Guidelines for Resuscitation 2010 Section 4. Adult advanced life support, Resuscitation 81, 1305-52
2. UK Resuscitation Council, Advanced Life Support Guide. 5th Ed., London: UK Resuscitation Council 2010
3. P Michalek, W Donaldson, L Theiler. The use of i-gel® in anaesthesia - Facts and fiction in 2013. Trends in Anaesthesia and Critical Care 2013 Oct; 3(5):246-251
4. L Theiler, M Gutzmann, M Kleine-Brueggene, N Urwyler, B Kaempfen, R Greif. i-gel® supraglottic airway in clinical practice: a prospective observational multicentre study. British Journal of Anaesthesia 2012 Dec; 109(6):990-5
5. M Kleine-Brueggene, L Theiler, N Urwyler, A Vogt, R Greif. Randomised trial comparing the i-gel® and Magill tracheal tube with the single-use ILMA® and ILMA® tracheal tube for fibre optic guided intubation in anaesthetised patients with a difficult airway. British Journal of Anaesthesia 2011 Aug; 107(2):251-7
6. D Haske, B Schempf, G Gaier, C Niederberger. Performance of the i-gel® during pre-hospital cardiopulmonary resuscitation. Resuscitation 2013 Sep; 84(9):1229-32
7. RM Beringer, F Kelly, TM Cook, J Nolan, R Hardy, T Simpson, MC White. A cohort evaluation of the paediatric i-gel® airway during anaesthesia in 120 children. Anaesthesia 2011 Dec; 66(12):1121-6
8. DA Gabbott, R Beringer. The i-gel® supraglottic airway: A potential role for resuscitation? Resuscitation 2007; 73(1): 161-2
9. P Michalek and W Donaldson (Edited by). The i-gel® supraglottic airway. Nova Science Publishers, 2013
10. RM Levitan, WC Kinkle. Initial anatomic investigations of the i-gel® airway: a novel supraglottic airway without inflatable cuff. Anaesthesia 2005; 60(10):1022-1026



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK  
 T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



The manufacturer Intersurgical Ltd is certified to ISO 9001:2015, ISO 13485:2016 and ISO 14001:2015

**Please think before you print**  
 Save energy and paper.  
 If you must print this information sheet please print it double sided.

IS6.3 • Issue 14 03.20